

GRIP™
by Acrostak



Full control over slippery and calcified lesions


Optimal for the treatment of **intra-stent restenosis**, **heavily calcified lesions** and for **stent post-dilation**.

Unique Knobs Design

- 4 lines of knobs on the surface of the balloon anchor the balloon to the lesion during balloon inflation
 - ▶ **Reduced risk of dissection** at the stent margins of ISR, and in delicate areas such as bifurcations, ostial lesions, and short calcified lesions
- 4 lines of knobs designed to crush the hard plaque remaining between the struts of the stent
 - ▶ **Optimal plaque crushing** for superior post-dilation

Outstanding Material Strength

- High pressure balloon with **RBP of 22 atm**
- **Ultra-resistant to puncture**
- **Truly tapered tip** for optimal crossability
 - ▶ **Securely reach and dilate** heavily calcified lesions without compromising **crossability**

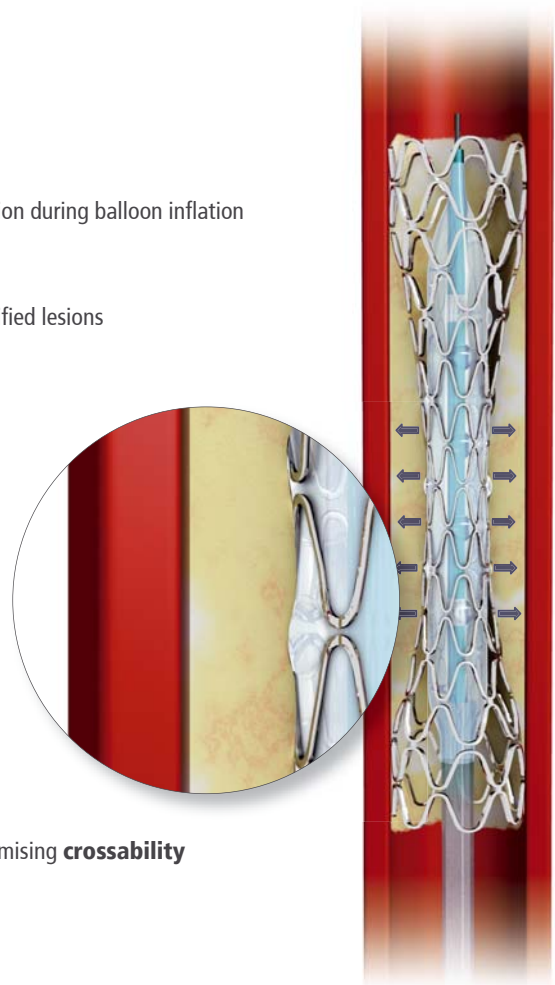


CE 0481

© ACROSTAK | All rights reserved.

Acrostak
Devices for interventional applications

GRIP™

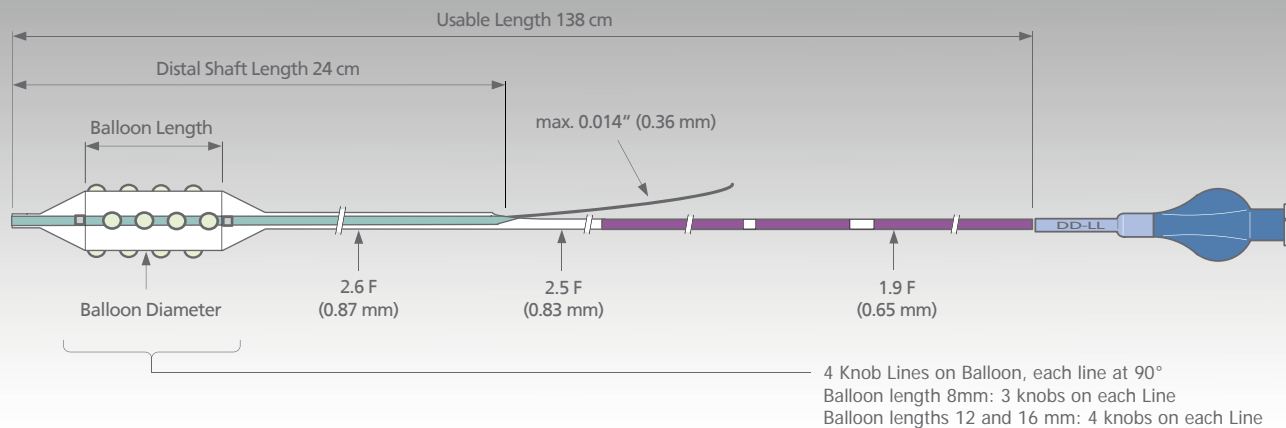

COMPLIANCE CHART

Pressure (atm)	Ø2.5 10/15/20	Ø3.0 10/15/20	Ø3.5 10/15/20	Ø4.0 10/15/20
6	2.30	2.76	3.19	3.65
8	2.37	2.84	3.31	3.78
10	2.45	2.94	3.44	3.91
11 ^N (1115 kPa)	2.50	3.00	3.50	4.00
12	2.54	3.05	3.56	4.07
14	2.62	3.14	3.67	4.19
16	2.68	3.22	3.76	4.29
18	2.73	3.27	3.83	4.37
20	2.77	3.33	3.90	4.45
22 [*] (2229 kPa)	2.81	3.38	3.97	4.52

* Rated Burst Pressure - do not exceed

^N Nominal

ORDERING INFORMATION

Reference no.	Diameter (mm)	Length (mm)
GRIP™		
250080340	2.5	8
250120340	2.5	12
250160340	2.5	16
300080340	3.0	8
300120340	3.0	12
300160340	3.0	16
350080340	3.5	8
350120340	3.5	12
350160340	3.5	16
400080340	4.0	8
400120340	4.0	12
400160340	4.0	16

Not authorised for sale in the U.S.



Worldwide distribution by:

Acrostak International
 Route de Malagnou 6
 CH-1206 Geneva
 SWITZERLAND

sales@acrostak.com
 Phone +41 22 347 68 70
 Fax +41 22 347 68 71

Manufactured by:

Acrostak (Schweiz) AG
 Stegackerstrasse 14
 CH-8409 Winterthur
 SWITZERLAND

info@acrostak.com
 Phone +41 52 233 95 51
 Fax +41 52 233 95 50

www.acrostak.com

Distributed by:

